



# FitMax™ CPAP Total Face Mask

for use with CPAP and bilevel NIV treatment




Critical Care ▪ FitMax™ CPAP Total Face Mask



### FitMax™ CPAP Total Face Mask

The FitMax CPAP total face mask is an adult (>30kg) single patient use NIV patient interface for use with CPAP and bilevel treatment within the hospital environment. The design of the mask helps to eliminate the nasal bridge issues that can be associated with traditional oronasal masks. The FitMax is easy to fit and is available in 3 sizes for the majority of patients.

Two colour coded elbow options are provided for different clinical situations – an anti-asphyxia elbow with leak port for use with single limb breathing systems and a standard elbow for use with dual limb breathing systems.

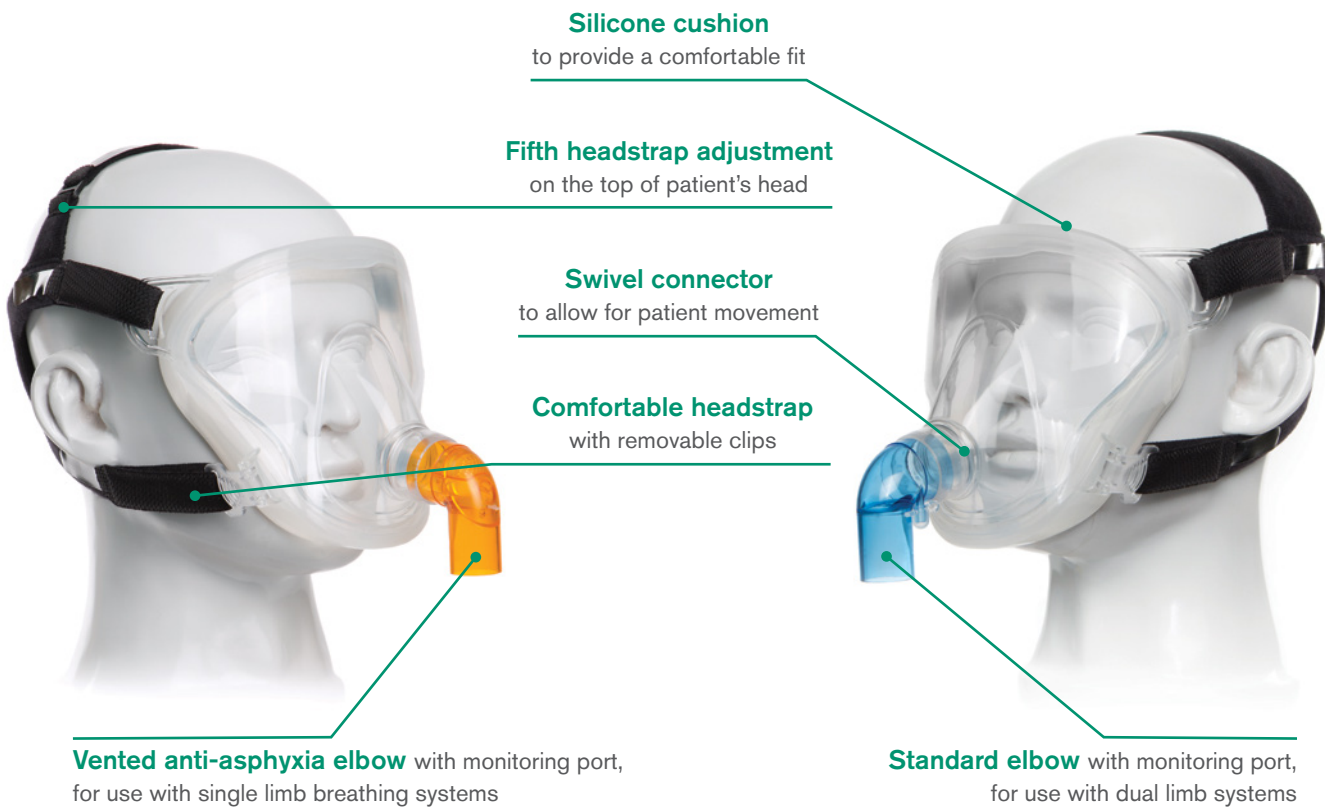

Code	Description	Colour coded elbow	Box Qty.
<b>For use with single limb breathing systems</b>			
2290020	FitMax CPAP total face mask, with vented anti-asphyxia elbow, head strap, medium	●	10
2290030	FitMax CPAP total face mask, with vented anti-asphyxia elbow, head strap, large	●	10
2290040	FitMax CPAP total face mask, with vented anti-asphyxia elbow, head strap, extra large	●	5
<b>For use with dual limb breathing systems</b>			
2290070	FitMax CPAP total face mask, with standard elbow, head strap, medium	●	10
2290080	FitMax CPAP total face mask, with standard elbow, head strap, large	●	10
2290090	FitMax CPAP total face mask, with standard elbow, head strap, extra large	●	5

Make an enquiry

The FitMax masks are manufactured by Besmed Health Business Corporation, Taiwan and distributed by Intersurgical Ltd



Intersurgical Benelux B.V. Vluchtoord 28, 5406 XP, Uden, Nederland  
T: +31 (0)413 243 860 F: +31 (0)413 243 869 info@intersurgical.be be-fr.intersurgical.com



Le fabricant Intersurgical Ltd est certifié ISO 9001:2015, ISO 13485:2016 et ISO 14001:2015

**Pensez-y avant d'imprimer**  
Economisons les ressources.  
Si vous devez imprimer ce document, faites-le recto verso.

IS1.52 BEFR • Issue 1 01.21